


# 緊急斷離程序 Emergency Disconnect Procedure

NxSTAGE  
供病人使用  
導管

1

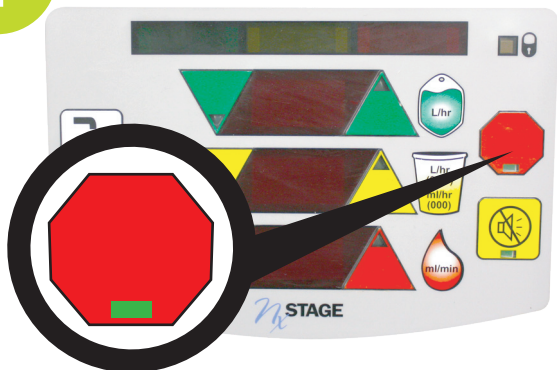


按血泵的按鈕  
血泵將停動，小長方格將轉成  
全綠色。

2




使用在急救袋中  
找到的用品...

3

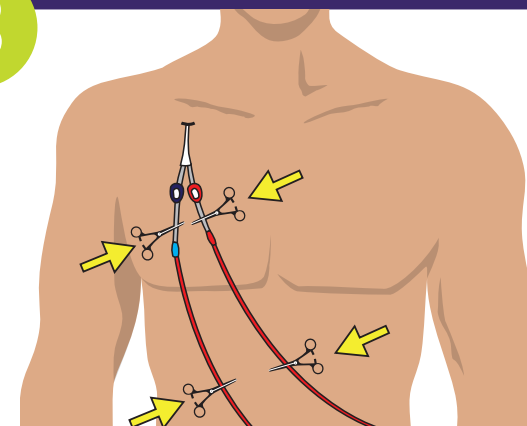

關好導管和血路管上的夾子。  
在連接位下方和上方加放  
安全夾。

4

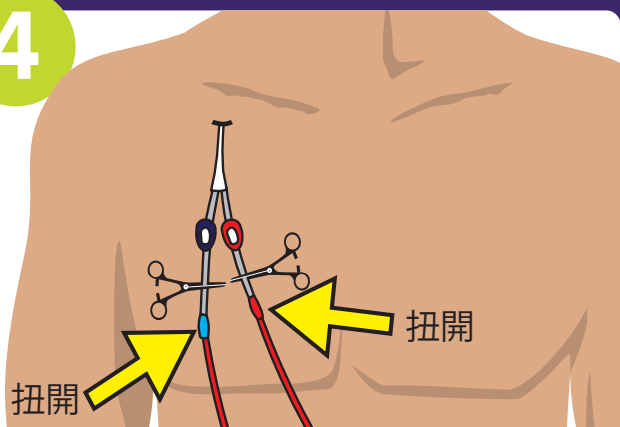

扭開  
扭開  
扭開連接位，  
將導管和血路管斷離。

5



攜同急救袋  
前往安全的地方

Produced by **BCRenal**  
Provincial Health Services Authority

**fraserhealth**  
Better health. Best in health care.

**Interior Health**

**northern health**  
the northern way of caring

**Providence Health Care**

**Vancouver Coastal Health**  
Promoting wellness. Ensuring care.

**island health**

Provincial Health Services Authority

**BRITISH COLUMBIA**